


"AMC-series" solves your problems on preparations of Cross-Sectional SEM specimens. No need to waste your time!

◆ Simple operations
Adjust→Teach→Scribe (Auto.)→Split
The target is skipped and never damaged.

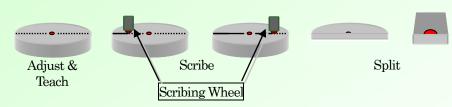

- ◆ CPU aided control of scribing (JPN Patent no. 3361024)

 AMC scribes both flat and curved specimens.
- ◆ e. g. AMC can split glass (9 mm thick) smooth.
 With pliers ABK-5: Aluminum Arm Type
 (for specimens > 5 mm thick)




AMC-7600

For Preparations of Cross-Sectional SEM Specimens.



AMC-7800



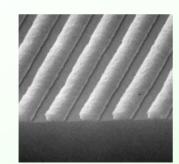

Cross-Sectional SEM image of mask glass



Specimens of glass

	AMC-7600	AMC-7800	AMC-8000
Specimens	Wafers, Photomasks, TFT-LCDs, Quartz, SiC, Sapphire, etc.		
Max. size	7"□ (8"φ)		
Thickness	~ 0.6 ~ 10 mm		
Scribing accuracy	± 15 μm		± 2.5 μm
Stage Movement	$x = 50 \text{ mm}, y = 230 \text{ mm}, \theta = 360 \text{ deg}.$		
Microscope (Magnification)	Stereoscope with a zoom lens $(\times 15 \sim \times 90)$ (High Mag. Type: $\times 45 \sim \times 280$)	CCD camera with a zoom lens $(\times 25 \sim \times 150 \text{ on a } 19 \text{ inch monitor})$	CCD camera with a zoom lens (×90~×580 on a 19 inch monitor)
Illuminator	Halogen lamp with a double light guide	Coaxial consisting of white LEDs	
Size	$W445 \times D800 \times H550 \mathrm{mm}$	$W445 \times D535 \times H620 \text{ mm}$	$W550 \times D650 \times H570 \mathrm{mm}$
Weight	~ 53 kg	~ 50 kg (without a LCD)	~ 78 kg (without a LCD)
Requirement	AC 100 V, ≤200W, 50/ 60 Hz		

Warranty: 1 year Custom Specifications are available by request.

Specifications are subject to change without notice.

AMC-series Options







A kind of Pliers to split wafers and glass (\leq 5 mm thick) easily

• Pliers ABK-5: Aluminum Arm Type

For specimens (> 5 mm thick)

- LED lamp with double light guides (AMC-7800/-8000)
- Tungsten Carbide Scribing Wheel

(a vertical angle of 105 deg./ 120 deg./ 135deg.)

• Diamond Scribing Wheel

(a vertical angle of 105 deg./ 120 deg./ 135deg.)




522-10, Kitano-machi, Hachioji-shi, Tokyo 192-0906, Japan Tel.: +81-42-646-4466(Rep.) Fax.: +81-42-656-1404 URL: http://www.apco-jp.com

E-mail: info@apco-jp.com

For Preparations of Cross-Sectional SEM Specimens

Precision Scriber and Pliers AMC-series

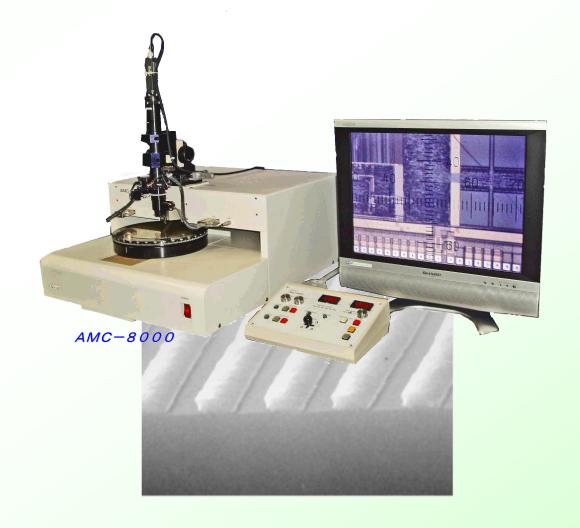

APCO Ltd.